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
### RICHWAY ™BIOMAT<sub>MX</sub>

Unlike an electric blanket, the Biomat does not get hot when it is turned on. The heat comes from friction generated within the body through the effects of Far Infrared Rays. It warms the inside of your body and the part of the Biomat making contact with your skin gets warm.

Should the temperature of any part of the Biomat exceed the set temperature, the sheet current is cut off by thermal sensors distributed throughout the Biomat. When the Biomat reaches the selected temperature, then the Negative Ions come on and will stay on until the Biomat needs to raise the temperature. It comes with an external control box which allows you to adjust the desired temperature and running time (between 2, 4, 8, or 12 hours). As an added safety feature, when the controller is set to 122°F or higher, it will automatically step down to 113°F after 4 hours.

## 6.5 Micron Infrared Ray

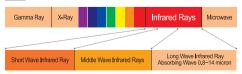Far Infrared Rays penetrate 6 inches into the body. Therefore, these rays not only benefit the muscles on the surface of the body, but all cells including blood vessels, lymph glands, and nerves in the deepest parts of our body.

## ™BIOMAT<sub>MX</sub> SPECIFICATIONS

- 1 The Super Fiber materials generate Far Infrared Rays (6-12 microns) that the human body can easily absorb.
- A special layer of double fiberglass and carbon completely cuts off harmful

## What does the Biomat EMF interceptor do?

The Biomat EMF Interceptor acts in the same way as a lightning rod, safely conducting the

EMF energy away to ground. With this feature.

you can effectively block exposure to EMF for yourself and your family.

BIOMAT EMF, EMC TEST (Tested by FCC L.)

General Electric device Intercepting System

fall well within safety levels established by the FCC (United States) VED (German) and CISPR (World)

**Biomat EMF Blocking Levels** 

# RICHWAY BIOMAT

**Licenses and Certificates** 

U.S. FDA

U.S. Food and Drug Administration(FDA) edical device (510K 072534)

The Ministry of Food and Drug Safety(Korea) approved medical device QUEEN Item Permit No. 13-1017 / SINGLE Item Permit No. 13-1018 PRO Item Permit No. 13-279 / MINI Item Permit No. 13-234

### J.ELECTRIC SAFETY

itted by the Japane Safety Authority (09-09087182)

Good Manufacturing Practice (GMP) approved manufacturer (No. 3965)

#### ISO 13485

Approved as an international medical device manufacturer (08 0037 SJ)

### ISO 9001

Approved as an international electrical equipment manufacturer (SKN1454/9474D)


International Electrical Safety Approval (Eui-si-2009-0555)

Korean Electrical Safety Corporation (KESCO)

PRO Item Permit No. HH071485-13002A MINI Item Permit No. HH071485-13001A

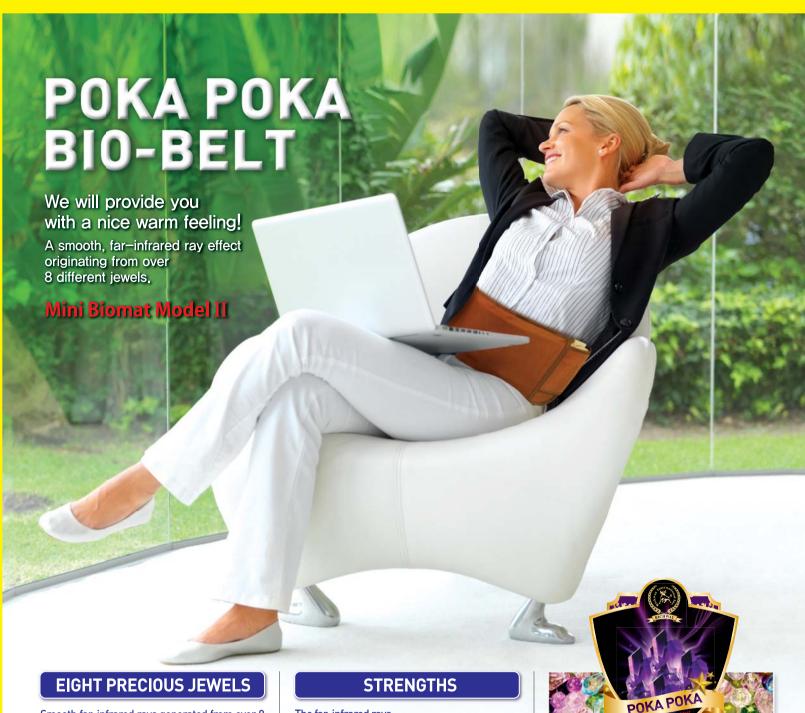
Smooth far-infrared rays generated from over 8 different types of jewels including Amethyst, Tourmaline, Green Jade, Crystal, Citrin, Topaz, Tiger Eye, and Elvan

- ► AMETHYST Known as the jewel of healing.
- ▶ GREEN JADE This jewel generates deep, smooth far-infrared rays.
- ► TOURMALINE The tourmaline generates natural anions.
- ► CRYSTAL This jewel contains strong in-body purification functions.
- ► CITRIN This jewel stabilizes the mind.
- ► TOPAZ This jewel generates strong far-infrared rays.
- ► ELVAN This jewel discharges bodily waste through the deep penetrating effect of far-infrared rays.

The far-infrared rays generated from over eight different jewels smoothly permeates the body.

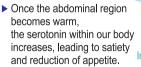- It is very convenient, as ailing body parts such as the back, waist, and knees can be treated easily without need of assistance.
- ► A pleasant and refreshing feeling comes from using the Poka Poka Bio-Belt!

  The belt can be used according to individual.

The belt can be used according to individual lifestyles-while relaxing and listening to your favorite music, watching DVD movies, working on the computer, and much more!

MEDIUM: 55.12"x5.12" / XL: 66.93"x5.12"

INFRARED RAY HEATING PART: 18.11" x 8.27"

• Because it is an electrical product, do not disassemble parts such as the controller and the belt, • Do not touch with wet hands, • It may cause electrical shock, • Restrict children or critical patients from using the product alone, • Make sure others supervise usage, • Avoid using at high temperature for pregnant women, • Handle with care. Do not throw or kick the product, • Never wash the product with cleansers, • Warming is set to be cut out after 1 hour at 60°C, It can be set repeatedly, • Because the controller needs ventilation, do not wrap it with vinyl or a blanket, • Do not place the belt in direct contact with skin. Use a towel, • Consult a doctor if you use a pacemaker.

\*There are individual differences in effectiveness.

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